

EXPERIMENT: To determine the amount of NaOH and Na_2CO_3 present in the same solution.

THEORY: When a mixture of NaOH and Na_2CO_3 is titrated against a standard HCl solution, colour of solution changes from yellow to pink (using methyl orange as an indicator) due to complete neutralization of both the alkalis at $\text{pH} \approx 4$. However, when the mixture is titrated using phenolphthalein as an indicator, the colour of the solution changes from pink to colorless due to complete neutralization of NaOH and half neutralization of Na_2CO_3 (i.e. upto the conversion of Na_2CO_3 to NaHCO_3) at $\text{pH} \approx 8$. The difference of two titre values gives the amount of HCl required for half neutralization of Na_2CO_3 while the diff. of first titre value and twice the second titre value gives the amount of HCl required for NaOH neutralization.

PROCEDURE:

(i) Standardization of HCl

1. Transfer 10ml of standard 0.1N Na_2CO_3 solⁿ in a clean conical flask using a pipette.
2. Add 2 drops of methyl orange indicator. Titrate the solution against HCl from the burette.
3. The colour of the solⁿ changes from yellow to pink (end point).
4. Note the volume of the solution used and repeat the titration at least 4 times and take the mean of closely related readings (V_0).

Khushi.

Teacher's Signature : _____

Expt. No. 1


Date: 29/09/20

EXPERIMENT: To determine the amount of NaOH and Na_2CO_3 present in the same solution.

APPARATUS: Pipette, burette, beakers, conical flask, funnel, burette stand and clamp.

CHEMICALS: Sodium carbonate (Na_2CO_3), sodium hydroxide (NaOH), hydrochloric acid (HCl), methyl orange and phenolphthalein.

CHEMICAL REACTIONS:

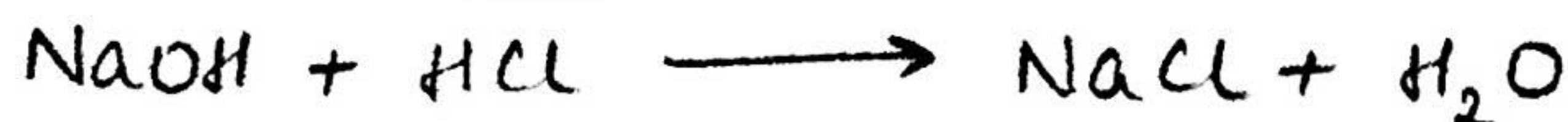

CHEMICAL STRUCTURES:

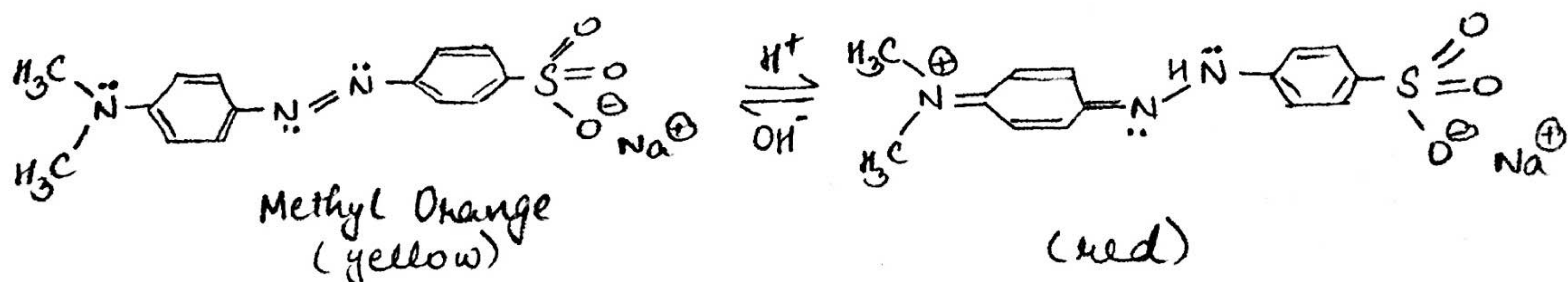


Fig. 1

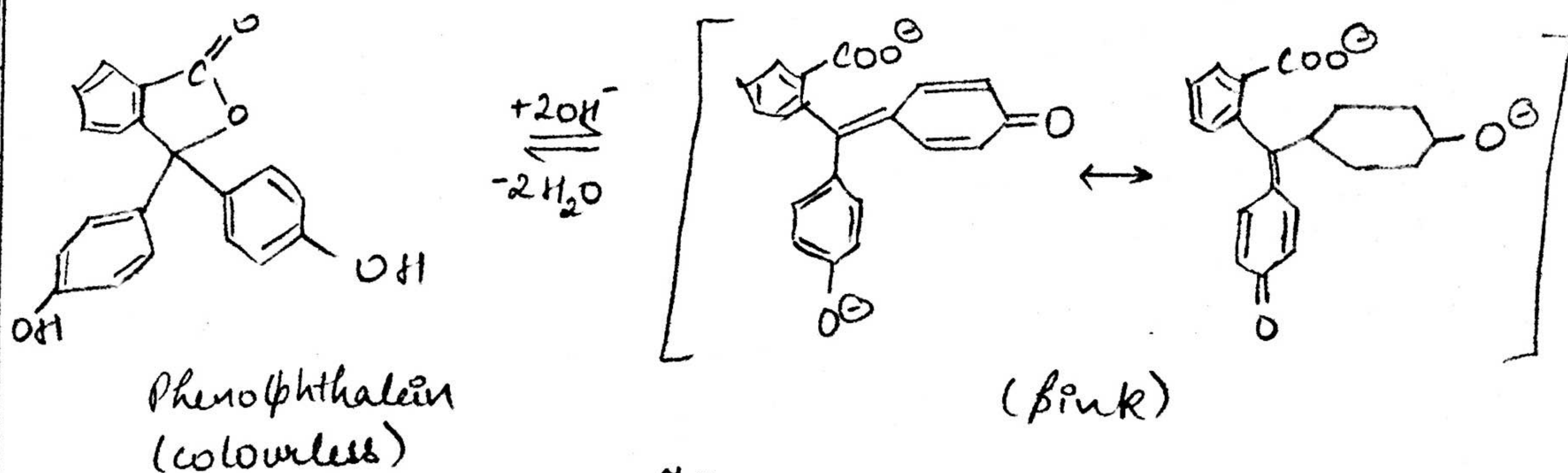


Fig. 2

INDICATOR: Methyl Orange and Phenolphthalein

Khusli

ii) Determination of NaOH and Na_2CO_3 content

1. Transfer 10ml of mixture of alkali solⁿ into a conical flask.
2. Add 2-3 drops of phenolphthalein indicator. The solution becomes pink in colour.
3. Note the initial reading of HCl from the burette (A). Titrate the solⁿ with standard HCl while the solⁿ becomes colorless (B).
4. Note the titre value and this is the phenolphthalein end point P. To the same solⁿ, add 2-3 drops of methyl orange indicator and continue the titration with HCl, until a sharp colour change occurs from yellow to red at the end point (C).
5. This titre value i.e., the total volume of HCl run down from the beginning of the experiment to the methyl orange end point is noted and this is the methyl orange end point M.

GENERAL CALCULATIONS:

Standardisation of HCl

Volume of alkali solution (Na_2CO_3) taken = 10ml

Normality of $\text{Na}_2\text{CO}_3 = N_2$; Volume of HCl used = V_0

Using the normality equation -

$$N_1 \times V_0 = N_2 \times 10$$

$$N_1 = N_2 \times (10/V_0)$$

Determination of Na_2CO_3 and NaOH

ii) Determination of NaOH

Equivalent weight of NaOH = 40

Hence 1L of 1N HCl = 40g of NaOH

Normality of HCl used = N_1

V_2 ml of N_1 HCl = $40 \times (V_2/1000) \times N_1 = y$ gm of NaOH

Khushti

Teacher's Signature: _____

OBSERVATIONS:

(i) Standardization of HCl solution

Volume of 0.1N Na_2CO_3 solution taken for each titration = 10ml

Sn. No.	Burette Reading (mL)		Volume of HCl used (mL)
	Initial	Final	
1.	0.1	10.0	9.9
2.	0.1	10.0	9.9
3.	0.1	10.0	9.9

Mean Value of HCl used (V_0) = 9.9 mL

(ii) Determination of NaOH and Na_2CO_3 in the mixture

Volume of mixture of NaOH and Na_2CO_3 solution taken for each titration = 10ml

Sn. No.	Burette Reading (mL)			Volume of HCl used (mL)	
	Initial Reading (A)	Colorless with phenolphthalein (B)	Reddish color with methyl orange (C)	$P = B - A$	$M = C - A$
1.	0.1	14.0	16.5	13.9	16.4
2.	0.1	14.0	16.5	13.9	16.4
3.	0.1	14.0	16.5	13.9	16.4

Mean volume of HCl used for P = 13.9 mL

Mean volume of HCl used for M = 16.4 mL

As P corresponds to $\frac{1}{2}$ neutralization of Na_2CO_3 and complete neutralization of NaOH, \therefore Half of $\text{Na}_2\text{CO}_3 = M - P = 16.4 - 13.9 = 2.5 \text{ mL}$

So, Volume of HCl required for neutralization of $\text{Na}_2\text{CO}_3 = 2(M - P) = V_1 = 5 \text{ mL}$

Volume of HCl required for neutralization of NaOH = $M - 2(M - P) = 2P - M = 11.4 \text{ mL}$

Khushi.

This is the amount of NaOH present in 10ml of the given alkali mixture solution.

$$\text{Strength of NaOH} = y_1 \times (1000/10) = 100 \times y_1 \text{ gm/L.}$$

(ii) Determination of Na_2CO_3

Equivalent weight of $\text{Na}_2\text{CO}_3 = 53$

Hence 1L of 1N HCl = 53 gm of Na_2CO_3

Normality of HCl used = N_1

$$V_2 \text{ ml of } N_2 \text{ HCl} = 53 \times (V_1 \times 1000) \times N_1 = y_2 \text{ gm of } \text{Na}_2\text{CO}_3.$$

This is the amount of Na_2CO_3 present in 10ml of the given alkali mixture solution.

$$\text{Strength of } \text{Na}_2\text{CO}_3 = y_2 \times (1000/10) = 100 \times y_2 \text{ gm/L.}$$

RESULTS: The given alkali mixture contains NaOH = $100 \times y_1 \text{ gm/L}$
= 4.6 gm/L

The given alkali mixture contains $\text{Na}_2\text{CO}_3 = 100 \times y_2 \text{ gm/L}$
= 3 gm/L

PRECAUTIONS: 1. Rinse the pipette with the solution to be transferred to titration flask.

2. Rinse the burette with the solution to be taken/filled in the burette.

3. Upper meniscus to be read for colored solutions.

4. Lower meniscus to be read for colourless solutions.

Khusli.

Teacher's Signature: _____

CALCULATIONS:

Standardization of HCl

using normality equation,

$$N_1 \times V_1 = N_2 \times V_2$$

$$N_1 = \frac{0.1 \times 10}{9.9} = 0.101 \text{ N}$$

Determination of Na_2CO_3 and NaOH

(i) Determination of NaOH

$$y. \text{ gm of NaOH} = 40 \times \frac{11.4}{1000} \times 0.101 = 0.046 \text{ gm}$$

$$\text{Strength of NaOH} = 0.046 \times 100 = 4.6 \text{ gm/L}$$

(ii) Determination of Na_2CO_3

$$y. \text{ gm of Na}_2\text{CO}_3 = 53 \times \frac{5}{1000} \times 0.101 = 0.03 \text{ gm}$$

$$\text{Strength of Na}_2\text{CO}_3 = 0.03 \times 100 = 3 \text{ gm/L}$$

RESULT:

Amount of NaOH present in 1L of alkali mixture = 4.6 gm/L

Amount of Na_2CO_3 present in 1L of alkali mixture = 3 gm/L