

EXPERIMENT: To determine pK_a value of acetic acid by pH metric titration.

THEORY: A pH meter will be used to follow the titration of an unknown weak acid, $HA(aq.)$ with sodium hydroxide, $NaOH(aq.)$. The weak acid has a concentration around $0.1M$. The result of the pH versus volume of $NaOH$ plot is 'S' shaped curve which is not as steep as the one arising from the titration of a strong acid. The equivalence point (this time) will be at alkaline pH (not 7 as in strong acid vs strong base). From the equivalence point, the concentration of an unknown acid HA is found. In addition, the acid constant K_a can be determined.

PROCEDURE:

- Titration of unknown HA with standard $NaOH$
1. Calibrate the pH meter with the standard buffer solution of $pH = 4$ or 9 , then rinse the glass electrode and immerse it in the beaker. Position the burette so that the titrant can be easily added.
 2. Pipette out $50mL$ of acetic acid into a clean beaker, dip the glass electrode. Record the pH.
 3. Initially, add $0.5mL$ of $0.1NaOH$ solution at a time, record the pH (after each addition), until the pH change is more than $0.2 - 0.3$ units, then start adding $0.2mL$ of $NaOH$ each time (i.e., near to the equivalence point, decrease the volume of $NaOH$ added) so that the change in pH is small enough.

Khushi

Teacher's Signature: _____

Date: 24/11/25


Expt. No. 8

EXPERIMENT: To determine pK_a value of acetic acid by pH metric titration.

APPARATUS: Pipette, burette, beakers, funnel, burette stand, clamp, pH meter and glass electrode.

CHEMICALS: Sodium Hydroxide (NaOH) and acetic acid (CH_3COOH).

CHEMICAL REACTIONS:

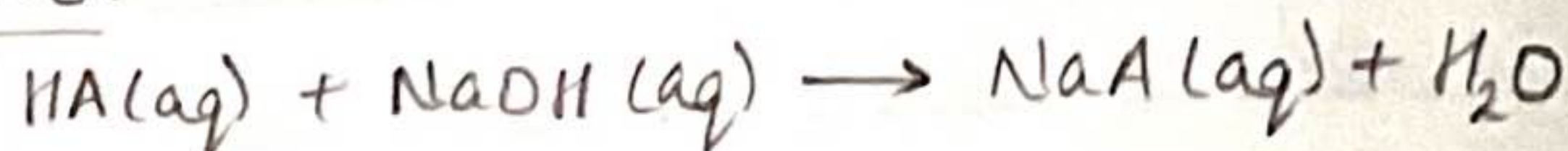

$$K_a = \frac{[H_3O^+][A^-]}{[HA]}$$

$$pH = pK_a + \log \frac{[salt\ form]}{[acid\ form]}$$

Henderson-Hasselbalch Equation

DIAGRAM:

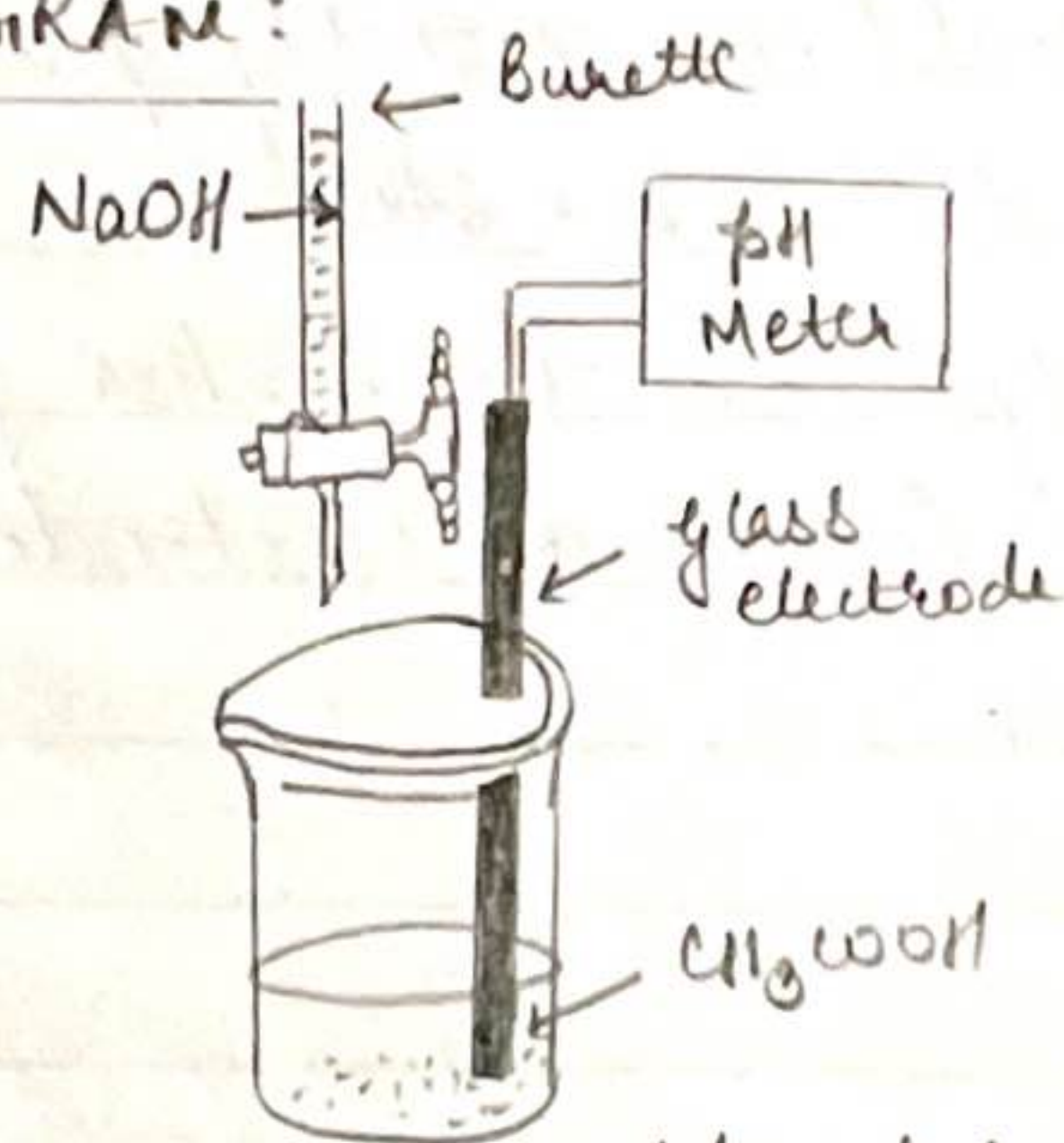


Fig. 1: Apparatus

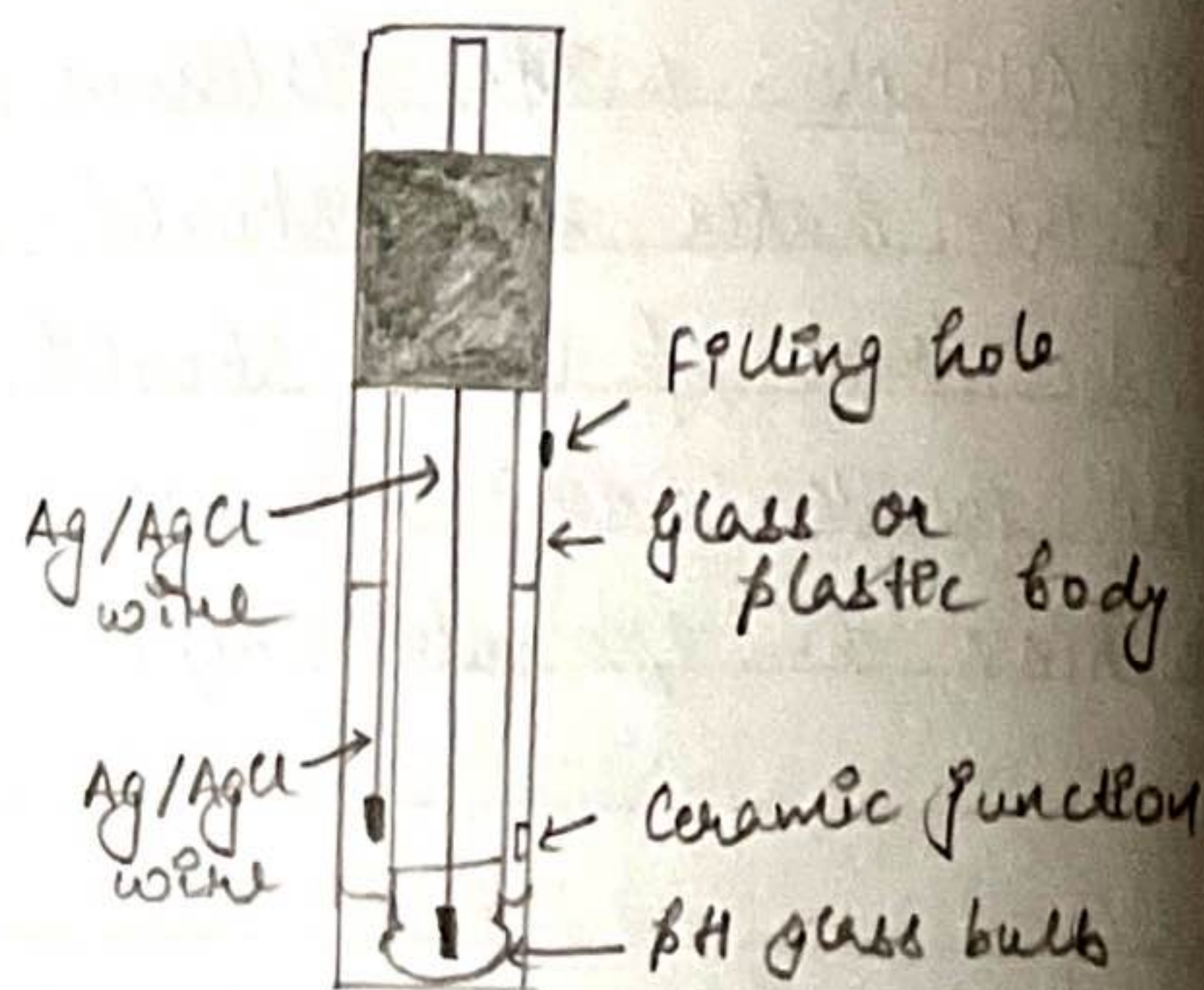


Fig. 2: Glass Electrode

to yield a good shape of plot.

4. After the rapid change in pH (after the equivalence point), the volume of NaOH may again be increased to 0.5 ml per addition. Make at least 10 more additions after the equivalence point so that the region with the plateau can be plotted.
5. pK_a is determined by explaining the titration curve. The negative log of K_a is pK_a and is same as the pH at half the volume of equivalence point.
($pH = pK_a$ when logarithmic term is zero which is zero once $[salt] = [acid]$. This is true at half equivalence point)

RESULT: The pK_a of acetic acid is 4.65.

- PRECAUTIONS:
1. Rinse the apparatus before experimenting
 2. The drops of NaOH should directly fall into the beaker from the burette and not on the wall of beaker or electrode.
 3. The glass electrode should be correctly placed in the beaker and should not be used as a stirrer.

Thushli

Teacher's Signature : _____


OBSERVATIONS:

Normality of standard $\text{NaOH} = 0.1\text{N}$

| S. No. | Volume of NaOH added from the burette (mL) | pH | S. No. | Volume of NaOH added from the burette (mL) | pH |
|--------|---|------|--------|---|-------|
| 1. | 0.5 | 2.84 | 16. | 8.0 | 5.45 |
| 2. | 1.0 | 3.76 | 17. | 8.5 | 5.61 |
| 3. | 1.5 | 3.95 | 18. | 9.0 | 5.96 |
| 4. | 2.0 | 4.09 | 19. | 9.5 | 9.35 |
| 5. | 2.5 | 4.23 | 20. | 10.0 | 12.00 |
| 6. | 3.0 | 4.35 | 21. | 10.5 | 12.35 |
| 7. | 3.5 | 4.43 | 22. | 11.0 | 12.51 |
| 8. | 4.0 | 4.55 | 23. | 11.5 | 12.63 |
| 9. | 4.5 | 4.64 | 24. | 12.0 | 12.69 |
| 10. | 5.0 | 4.74 | 25. | 12.5 | 12.75 |
| 11. | 5.5 | 4.84 | 26. | 13.0 | 12.81 |
| 12. | 6.0 | 4.91 | 27. | 13.5 | 12.86 |
| 13. | 6.5 | 5.04 | 28. | 14.0 | 12.89 |
| 14. | 7.0 | 5.17 | 29. | 14.5 | 12.93 |
| 15. | 7.5 | 5.30 | 30. | 15.0 | 12.94 |

Khushi.

GRAPH:

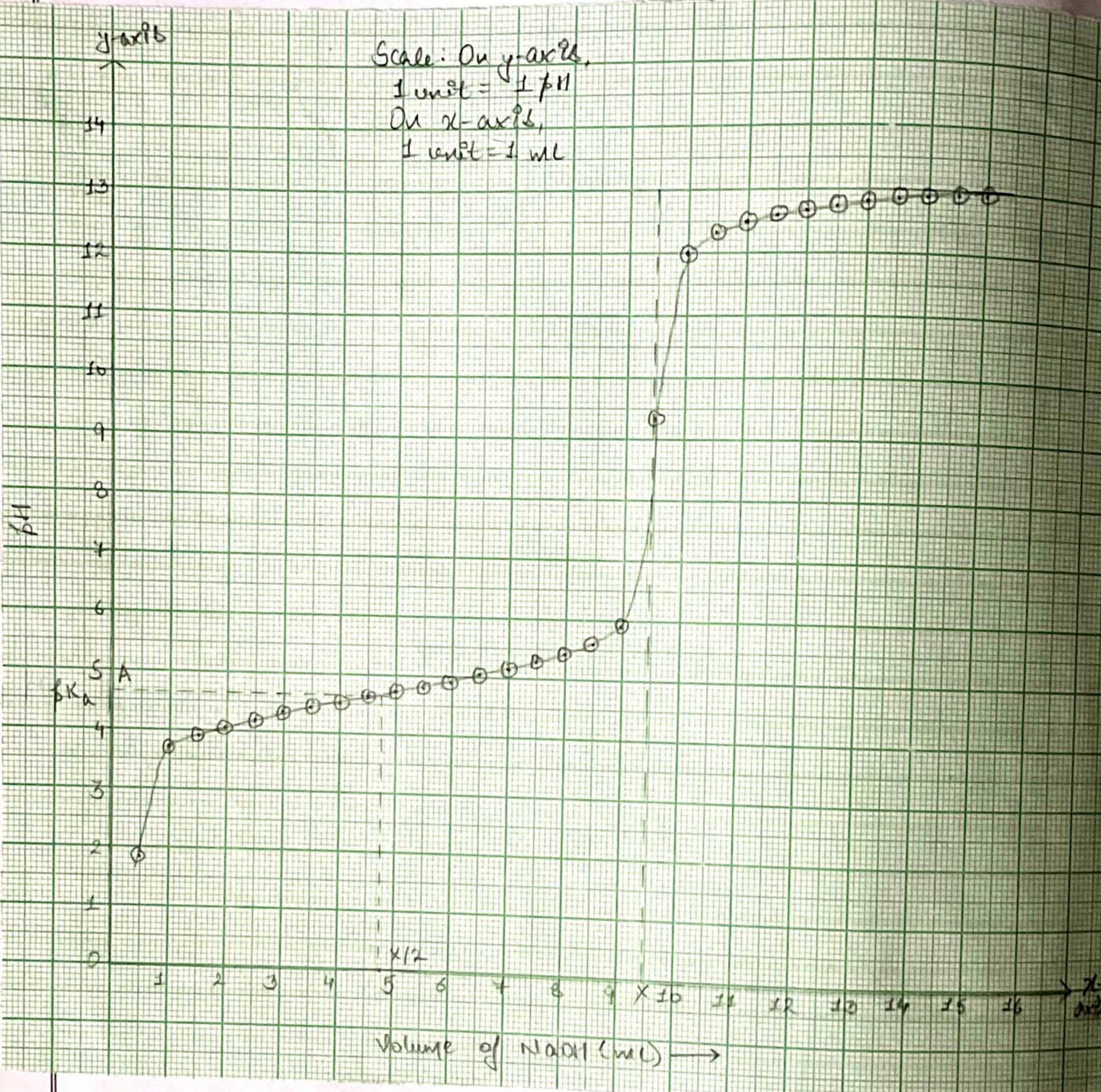

RESULT: The pK_a of acetic acid is 4.65.